
bruxane 2go

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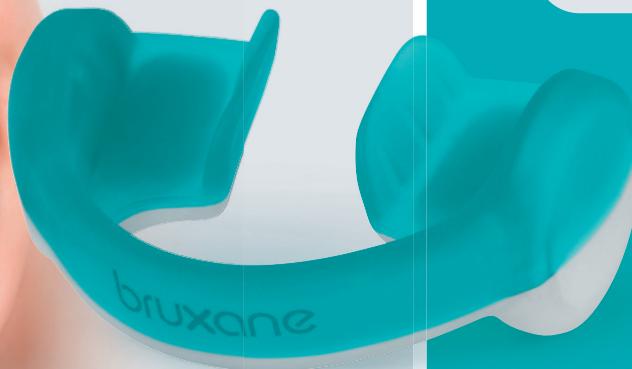Instruction manual

Mode d'emploi

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Fragebogen
innenliegend.

Questionnaire
inside.

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Inhalt

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Contents

Explanation of the Symbols _____	14
Explanation of the Signal Words _____	14
Intended Use _____	14
Safety Notices _____	15
Contraindications _____	16
Items Supplied _____	17
Product Overview _____	17
Before Initial Use _____	17
Use _____	18
Cleaning and Storage _____	18
Disposal _____	18
Product Data _____	19
Contact Details _____	19
Troubleshooting _____	20
FAQs and Further Information _____	20
Questionnaire to Monitor Your Personal Success _____	21

IMPORTANT NOTICES!

PLEASE KEEP FOR REFERENCE!

Before using the product for the first time, please read through these operating instructions carefully and keep them for future reference. They form an integral part of the product.

The manufacturer and importer do not accept any liability if the information in these operating instructions is not complied with.

Dear Customer,

We are delighted that you have chosen the **bruXane b2go**.

If you have any questions about the product, you can obtain further information and details of how to contact us on the website **www.bruXane.de** or by calling +49 (0) 6421 8079859*

You will also find these operating instructions on our homepage:

<https://www.bruXane.de/downloads>.

* Calls to German landlines are subject to your provider's charges

Explanation of the Symbols

The following symbols are used in the operating instructions, on the packaging, the product and any accessories:

Danger symbol: This symbol indicates dangers of injury. Read through the associated safety notices carefully and follow them.

Type BF (Body Floating)
In direct contact with patients for diagnosis, treatment or therapy. Must NOT have direct contact with the heart.

Manufacturer

Read operating instructions before use!

Electrical and electronic products must be properly disposed of in accordance with the WEEE Directive

Dispose of batteries properly.

Explanation of the Signal Words

WARNING

Warning of *possible* serious and life-threatening injuries.

NOTICE

Warning of damage to property.

Intended Use

The b2go dental splint aids the treatment of bruxism¹.

The b2go is designed, when there is corresponding pressure on the teeth, to produce a vibration which, by means of bone conduction to the ear, can also be perceived acoustically as a buzzing sound. With this biofeedback, the b2go subconsciously makes your body aware of its “grinding/clenching” malfunction so that it can immediately stop doing this.

- The b2go should only be worn by yourself for a short time until you have achieved the successful outcome you want.
- The maximum lifespan of the b2go is 8 weeks, with an average number of 100 grinding episodes per night.

¹ Bruxism = the subconscious grinding or clenching of teeth, usually at night.

- Do not wear the b2go for longer than 8 weeks. If you wear it for longer periods, this could cause elongation² of the teeth.
- The lifespan and durability of the material and the battery of the b2go may be shorter with very severe and frequent grinding or intentional clenching.
- The b2go is for personal use only and is not intended for commercial applications.
- The b2go should only be used as described in these operating instructions. Any other use is deemed to be improper.
- The warranty does not cover faults caused by incorrect handling, damage, assembly or attempted repairs. The same applies to normal wear and tear.
- Any dental splint can be bitten through if it is used excessively. No warranty can be offered against biting through it.

Safety Notices

WARNING – Danger of Injury

- **Danger of suffocation!** Keep children and animals away from the packaging material.
- b2go is **not** a toy! Keep it away from children and animals.
- **Danger of burns!** If the battery has leaked, do not allow your skin, eyes and mucous membranes to come into contact with the battery acid. If you touch battery acid, rinse the affected areas with abundant clear water straight away and seek medical assistance immediately.
- b2go must not be modified, manipulated or taken apart.

NOTICE – Risk of Material Damage

- Do not press any sharp objects into the b2go (pens, finger nails, etc.).
- Keep your b2go away from ultrasonic devices or microwave ovens.
- b2go can only vibrate when it is inside the oral cavity and the teeth are clenched against each other. If it vibrates constantly, stop using the b2go and go to the FAQ page or contact us by telephone.
- Consuming sleeping pills or sedatives or alcohol before wearing your b2go may reduce the effectiveness of any biofeedback and therefore cause increased grinding activity. This increased load may shorten the useful life of the material (biting through it) and the battery. No warranty can be offered against biting through it.
- Do not clean in the dishwasher.
- Do not use any denture cleansers, dips, toothpaste or other sharp or abrasive cleaners. These contain strong bleaching or cleaning agents which may attack the material.

² Elongation = minimal displacement of teeth, usually into the gap of the opposing jaw.

Contraindications

Read through the following list carefully. If you fall into one of the groups of people listed, do **not** wear the b2go or discuss wearing it with your doctor before you do so. If you are unsure, you should always consult your doctor.

- b2go must not be used by children (under 18 years of age).
- b2go may not be used by people with reduced physical, sensory or mental abilities or a lack of experience and / or knowledge.
- The following groups of people must not wear the b2go:
 - Pregnant and breastfeeding women
 - Asthmatics
 - Epileptics
- Do not wear the b2go if you have serious malocclusions and jaw abnormalities that make it impossible to wear a normal dental splint.
- Do not wear the b2go if you have large gaps in your teeth with at least two consecutive teeth missing in the area of the premolars.
- If your cheek area is too small (e.g. if a boiled sweet causes discomfort in the cheek area), you should refrain from wearing the b2go because you will find the splint uncomfortable and bothersome.
- If you have excessive sensitivity to touch (hyperaesthesia) in the cheek area, you should not wear the b2go.
- Do not wear the b2go if you know you are allergic to one of the ingredients of silicone.
- b2go must also not be used:
 - if you have serious difficulty breathing (e.g. asthma, colds, coughs)
 - while you are being treated with a defibrillator
 - if you are using flammable anaesthetics and gases
 - if you are under anaesthetic
 - in strong electromagnetic environments
- In the case of patients with the following conditions, b2go may only be used **following explicit prior consultation** with the doctor treating them:
 - in the case of body implants (excluding dental implants), for example heart pacemakers, defibrillators, hearing aids, insulin pumps
 - if there is an elevated risk of obstructive sleep apnoea syndrome
 - if there is an elevated risk of Upper Airway Resistance Syndrome (“UARS”)


Items Supplied

- bruXane b2go dental splint
- Storage box
- Operating instructions

Check the components for transport damage. If you find any damage, do not use the b2go but contact our customer service department.

Before using the b2go, remove all packaging materials and any labels or transport protection that may be present.

Product Overview

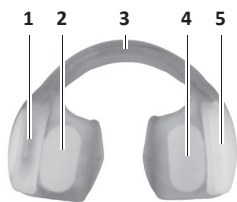- 1 Motor unit with indicator paper
- 2 Pressure sensor with indicator paper
- 3 Metal-reinforced link
- 4 Pressure sensor with indicator paper
- 5 Battery unit with indicator paper
- 6 Storage box

Before Initial Use

The b2go contains a flexible metal reinforcement in the anterior region (link between the two sensor areas) which allows you to adjust your b2go better to your own personal dentition. To do this, hold this link between your thumb and index finger close to the two sensor surfaces and bend the radius of the link narrower or wider - depending on whether you want a smaller or wider b2go shape.

Use

Do a visual inspection of your b2go **each time before you use it**:

- Each of the indicator papers must be white. A change in colour (red) may indicate a leak of the material or the battery.
- The silicone must not be damaged.

If either of these two cases arises, the b2go must not be inserted and must no longer be used.

1. Push the b2go onto your upper or lower jaw evenly using both hands.
2. Check that it fits correctly: The b2go must be positioned with the middle of the link at the mid-point of your front teeth.
3. Perform a functional test: Press your molars together (using roughly the force you would use to chew bread). This causes the b2go to vibrate. As soon as you reduce the pressure, the vibration stops. If the b2go does not vibrate, stop using your b2go and go to the FAQ page or contact us by telephone.
4. Wear your b2go while you are sleeping (even if you have a daytime nap), but for no more than 10 hours per day.
5. To extract the b2go, hold it on both sides and carefully take it out of your mouth.
6. Clean the b2go thoroughly after every use (see the “Cleaning and Storage” chapter).

Cleaning and Storage

Rinse the b2go under cold or lukewarm, clear water; use a little liquid soap if necessary. Then dry thoroughly with a lint-free cloth and place in the storage box. Always keep the storage box out of the reach of children and animals.

Disposal

Please dispose of the packaging material in an environmentally friendly manner and take it to a recycling centre.

The product is governed by the European Directive 2012/19/EU on waste electrical and electronic equipment (WEEE). Do not dispose of the product with built-in battery as normal domestic waste, but rather in an environmentally friendly manner via an officially approved waste disposal company.

If you have any questions on this, please contact your local waste disposal authority.

Product Data

Article number:	04212
Model number:	WB-01
Permissible ambient temperature:	Operation: 5 – 45 °C Storage and transport: 5 – 45 °C
Permissible humidity:	Operation: 30 – 80 % Storage and transport: 30 – 80 %
Permissible air pressure:	Operation: 500 – 1060 hPa Storage and transport: 500 – 1060 hPa
Dimensions:	approx. 2.3 x 2.4 x 8.5 cm (H x W x L)
Weight:	15 g
Electrical data:	Non-replaceable button cell CR2032, 3 V DC, 210 mAh
Product lifespan:	8 weeks
ID of operating instructions:	Z 04212 M bX V2 0720 uh

This product meets the requirements of the MDR Regulation EU 2017/745

Contact Details

Manufacturer within the meaning of the Medical Products Act and EU authorised representative according to Directive 93/42 EEC:

Hydas GmbH & Co. KG
 Hirzenhainer Straße 3, 60435 Frankfurt, Germany
 UDI DI: 4004128026337

Import:

bruXane GmbH
 Steinweg 43, 35037 Marburg, Germany
 Tel: +49 (0) 6421 8079859*

* Calls to German landlines are subject to your provider's charges.

Troubleshooting

Problem	Solution
Irritation of the oral cavity	Stop wearing your b2go until the irritation has subsided. Or try wearing the b2go on your other jaw. If the irritation persists, consult your dentist.
In the morning after wearing the b2go, you have the sensation “that your teeth do not fit together”	This sensation can occur in the morning, but normally disappears after a short time. If it persists, consult your dentist.
Allergic reaction	Stop wearing your b2go and consult your dentist.
Uncomfortable feeling while wearing / foreign body sensation	You should have got used to your b2go after a few nights. If you are using a dental splint for the first time, this may take a little longer.
Excessive salivation or dry mouth	You should have got used to your b2go after a few nights. If you are using a dental splint for the first time, this may take a little longer.
Waking up due to the biofeedback from the b2go	Patients who are very light sleepers may be woken at the moment of clenching or grinding. The body is then aware that you have been clenching/grinding. The body registers this and you will clench/grind less in the future. You should have got used to your b2go after a few nights.
The b2go falls out.	As the b2go is a standard size, it is normal for the dental splint to fall out of the mouth occasionally. However, the bruXane effect can have an impact after just a few hours of wearing it. The b2go can also be better adapted to your own jaw shape using the flexible link, as described in the “Before Initial Use” chapter.

FAQs and Further Information

You will find the latest questions and answers on our website: <https://www.bruXane.de/support/user/bruXane-2go>.

Questionnaire to Monitor Your Personal Success

To monitor the effectiveness of your b2go, we have developed the following questionnaire for you. Please rate the symptoms that apply to you before you first use your b2go and then after a few nights.

Questionnaire Before

Rating scale 0 to 10 (0 = not applicable, 10 = very strong)

Symptoms After Getting Up	Before
Headaches	
Tinnitus	
Impaired vision	
Aching facial muscles	
Tiredness in the jaw region	
Stiffness of movement of the jaw	
Rubbing sounds from the mandibular joint	
Clicking from the mandibular joint	
Aching mandibular joint	
Lower jaw blocked	
Shoulder pain	
Neck pain	
Back pain	
Fatigue / tiredness	
Other symptom that you ascribe to teeth grinding:	

Questionnaire After

Rating scale 0 to 10 (0 = not applicable, 10 = very strong)

Symptoms After Getting Up	Before
Headaches	
Tinnitus	
Impaired vision	
Aching facial muscles	
Tiredness in the jaw region	
Stiffness of movement of the jaw	
Rubbing sounds from the mandibular joint	
Clicking from the mandibular joint	
Aching mandibular joint	
Lower jaw blocked	
Shoulder pain	
Neck pain	
Back pain	
Fatigue / tiredness	
Other symptom that you ascribe to teeth grinding:	

My b2go helped me, what now?

Have you worn your b2go and noticed an improvement in your symptoms? Then **bruXane personal** is the right choice for you. To receive a bruXane personal that is custom-made for you and has a longer lifespan and is more comfortable to wear, please visit our bruXane website www.bruXane.de or ask your dentist.

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