Massage Wand

Item No. 56003

Instruction Manual

Please read all instructions carefully and retain for future use


(EN) Content

1. 1 x Massage Wand
2. 1 x Manual

Specifications

Dimensions	12.5 x 4.5 x 32cm
Material	ABS
Battery type	4 AA Batteries (Not included)

Instructions

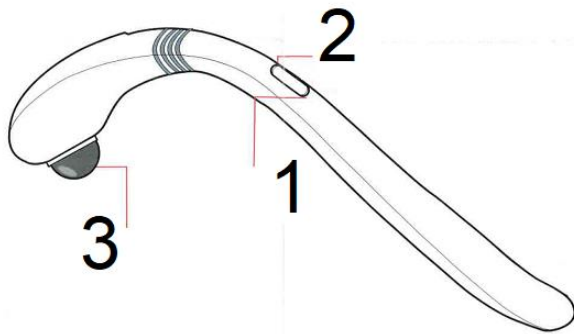1. On/Off
2. Velocity button
3. Massage node

Instructions (How to use the Massager wand)

1. Remove the unit from the box.
2. Remove any packaging from the product.
3. Remove the battery cover.
4. Insert 4 x AA Batteries (Not included) as per the polarity indicates.
5. Slide the switch to the Lower/High Vibration

Note: Do not leave product unattended when powered on. Remove the batteries when not in use.

Instructions (How to clean and maintain)

1. Only use a soft, dry cloth to clean.
2. Place the massager in a safe, cool and dry place when not in use.
3. After unplugging, clean the massager with a soft, damp cloth using NO abrasive cleaners.
4. Keep away from all solvents or harsh detergents.

Attention/Danger!:

Warning: Any user who is pregnant, has a pacemaker, suffers from diabetes, phlebitis, and/or thrombosis, heart problems or is at an increased risk of developing blood clots due to recent surgery, should consult their doctor before using a massaging device designed for home use.

Warning: To reduce the risk of fire or electric shock, do not expose this apparatus to rain or moisture. The apparatus shall not be exposed to dripping or splashing, and objects filled with liquids, such as vases, shall not be placed on apparatus.

Warning: Before use check the product for visible external damage. Don't use the product if it is damaged or has been dropped.


Precautions:

- Consult your doctor before using this product.
- Individuals with pacemakers should not use this item.
- This product may not be used on (small) children.
- This unit is not intended for use by persons with reduced physical, sensory or mental capabilities or lack of experience or knowledge. They should be given supervision and instruction in the use of the unit by a person responsible for their safety.
- Never use the product on open skin eruptions, swollen, or inflamed areas.
- Don't operate after taking painkillers, sedatives or alcohol. Consult your doctor if you are under these medications.
- Don't continuously operate for 15 minutes or longer. It is recommended to turn the product off and allow it to rest (cool) for 30 minutes before re-using. This will prolong the life of the product.
- The product is designed as a non-professional massager to soothe tired and arching muscles. Persist or worsen after use, contact your doctor immediately.

WEEE:

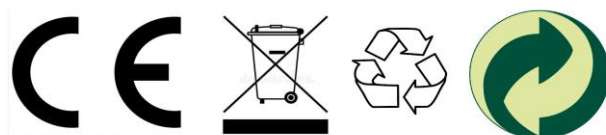This product bears the selective sorting symbol for Waste Electrical and Electronic Equipment (WEEE). This means that this product must be handled pursuant to European Directive 2012/19/EU in order to be recycled or dismantled to minimize its impact on the environment. The user has the choice to give this product to a competent recycling organization or to a retailer when buying a new electrical or electronic piece of equipment.

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