
FLOSSERPIK® POWER SPA

BEDIENUNGSANLEITUNG
INSTRUCTION MANUAL
MANUEL D'UTILISATION
MANUALE D'USO
MANUAL DE INSTRUÇÕES

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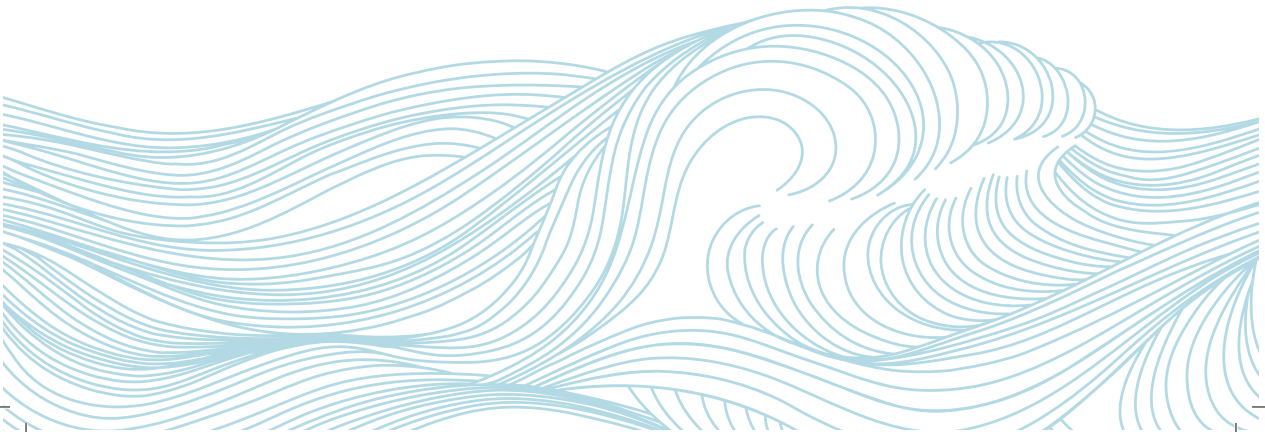For your safety and continued enjoyment of this product, always read the instruction book carefully before using and keep it in a safe place for future reference.
edel+white® and Flosserpik® are internationally registered trademarks of Scanderra GmbH, Basel, Switzerland. All rights and modifications reserved.

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DEUTSCH

CONGRATULATIONS

on your purchase of the edel+white® FlosserpiK® Power Spa from the Swiss Dental Expert brand, edel+white®.
With correct and consistent use, this water flossing system will keep your mouth fresh and help sustain your oral health indefinitely.

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FOR SAFETY REASONS

Warning signs instructions:

The following signs are used to classify and describe the degree of danger, injury and property damage caused by ignoring the sign description and improper use.

WARNING

Indicates a potential hazard that could result in serious injury or death.

Attention

Indicates a potential hazard that could result in minor injury or property damage.

The following marks are used to classify and describe the types of instructions to be followed.

This symbol is used to remind the user of the prohibited operation steps.

This symbol is used to remind users of the operation steps that must be followed in order to operate the product safely.

TIPS

This symbol is to improve operational steps.

Please note that your product may differ slightly in appearance to the illustrations on the following pages.

Safety Instructions:

Keep the product packaging out of reach of children.

Oral irrigators can be used by children aged 8 years and above and by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe manner and understand the hazards involved.

Children should be supervised to ensure that they do not play with the appliance.

Oral irrigators should be used only with water, or specific solutions as defined in the manufacturers instructions.

CAUTION:

Don't charge the product in humid surroundings.

Do not exceed the rating of the domestic power socket or wiring when charging. If the load of socket or wire exceeds the standard, it may cause fire.

If there are any abnormalities or failures, please stop using or charging the product to prevent electric shock, injury or fire.

 Abnormalities and failures include:
Broken device, deformed cable, abnormal heating, burning smell or sound during usage. In the event of any of these abnormalities, contact the customer service centre.

 When cleaning the power adaptor and charging cable, be sure to remove the power adaptor from the household power outlet.

 Do not repair the device and ensure proper maintenance to avoid accidents and injuries.

 Do not take the device apart

 This water flosser system should never be left unattended when plugged in.

 Do not place or store the unit where it can fall or be pulled into the tub, sink or toilet.

 Do not drop or knock the device

Never block the air openings of the jet stream

Keep the air openings free of lint, hair, and the like

Never drop or insert any object into any opening or into suction pipe

Do not remove the flexible tip on any of the attachments. If it separates for any reason, discard both the tip and shaft and replace them with new attachments.

Do not fill the tank with water warmer than 40°C

Fill reservoir with lukewarm water only unless otherwise specified by your dental professional.

 Do not overfill.

MEDICAL WARNING:

- If you have had recent oral surgery, consult your dentist before using the device.
 - When you start any new oral hygiene programme, your gums may bleed slightly. Consult your dentist if excessive bleeding continues after 10 days.
 - Do not use mouthwashes or rinses that contain iodine, hydrogen peroxide or tree oil.
 - Do not use while wearing oral jewellery. Remove prior to use.
 - Do not use if you have an open wound on your tongue or in your mouth.
 - If your physician or cardiologist has advised you to receive antibiotic premedication before a dental procedure, you should consult your dentist before using this oral hygiene aid.
 - Be sure that only the tip of the water jet touches your teeth, gums or tongue. Damage to those areas may occur if any other part of the handle touches them.
 - Do not direct water under the tongue, into the ear, nose or other delicate areas. This product can produce pressures that can damage sensitive areas.
-  See recommended technique for correct usage.
- Use only for oral hygiene cleaning. Do not use for any other purpose.
 - As with any new oral hygiene programme, individuals with heart problems should check with their physician prior to using new medical devices.

THE FLOSSERPIK® POWER SPA ORAL WELLNESS SYSTEM

Water flossing significantly contributes to a healthy mouth by helping to clean better all around the teeth and gum line. The FlosserpiK® Power Spa offers a double cleaning system, comprising a pulsating jet stream and toothbrushing with an UltraSoft Konex HD® brush.

WHY USE EDEL+WHITE® FLOSSERPIK®?

- PULSATING JET STREAM
- ULTRASOFT TOOTHBRUSHING

Daily irrigation helps reach places brushing simply cannot, promoting cleaner teeth and healthier gums. The result is an incredibly clean, fresh and healthy mouth. FlosserpiK® Power Spa is a powerful water pump, and the UltraSoft toothbrush helps remove debris around teeth and below the gum line. Water flossing with Power Spa reduces the risk of bacteria build-up in crevices and deeply cleans dental pockets and hard-to-reach spots, especially for people with implants, bridges and braces.

BEFORE USING IT FOR THE FIRST TIME

CHARGING


! For charging, make sure that the product and your hands are dry, then connect the USB charging cable to the power supply adaptor* and plug the power plug into the power supply socket (fig. 1).

Plug in the charging plug of the USB charging cable, and ensure plug is secure (fig. 2).

Red indicator is flashing after use or if the product can't be turned on. At this point, the battery is insufficiently charged and needs to be plugged in for charging (fig. 3).

- It takes 5-6 hours to fully charge.
- When charging, a green, slowly glowing indicator light is on.
- The green charging indicator light is on, once the charging is complete.
- When the battery is fully charged, the charging will be turned off to protect the battery from overcharging.


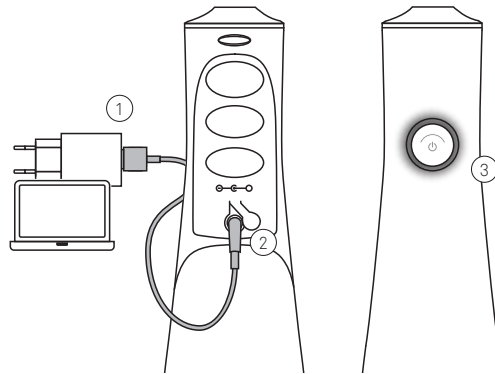
*adaptor not included

GETTING STARTED

ATTACHMENTS

- 1.) Jet attachment: suitable for rinsing away food residue.
- 2.) Tongue Scraper: suitable for cleaning the tongue.
- 3.) Subgingival attachment: recommended for people with periodontitis and other gum problems.
- 4.) Brush tip with UltraSoft Konex HD® bristles (Orthodontic attachment): standard tip recommended for people wearing braces or dentures.

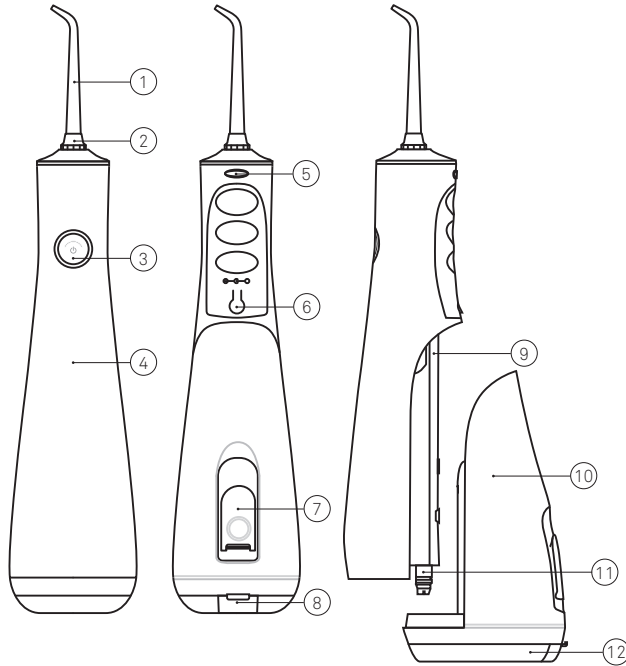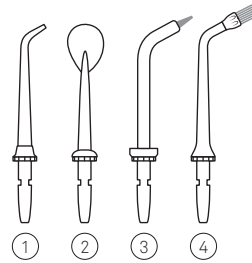- 1.) Attachment
- 2.) Attachment ring
- 3.) Power and selection button
- 4.) Main body
- 5.) Attachment release button
- 6.) Socket for DC charging
- 7.) Tank cover
- 8.) Hatch release button
- 9.) Suction pipe
- 10.) Water tank
- 11.) Filter
- 12.) Hatch cover

INSERTING AN ATTACHMENT

Insert the attachment into the centre of the handle (fig. 1). Gently try to pull the tip out to confirm the tip is locked into position (fig. 2).

REMOVING AN ATTACHMENT


Press the attachment release button (fig. 3) and remove the attachment at the same time (fig. 4).

FILLING THE TANK

Fill the water tank by opening the tank cover, while holding the device horizontally. Then close the cover of the tank (fig. 5).

 Fill water tank with liquid no warmer than 40°C.
Do not pour in salt or toothpaste.


 Wash the tank with water before first use.
Fill the tank with water, point the attachment towards the sink and press the power button to clean the device. Empty the tank fully to ensure the device is fully cleaned.

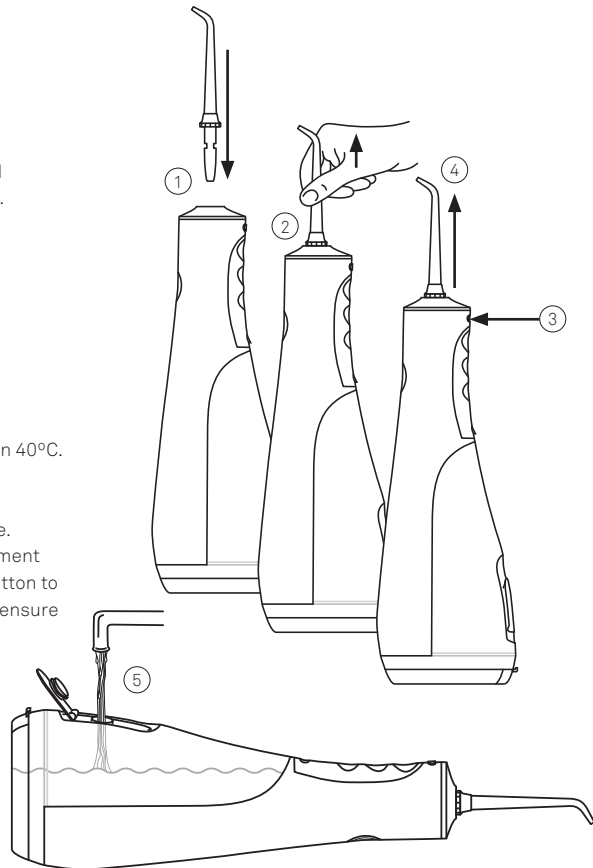HOW TO USE YOUR FLOSSERPIK® POWER SPA ORAL WELLNESS SYSTEM

⊘ Do not power on the device before placing the tip in your mouth.

1. Bend over the sink and hold the handle.
Put the attachment into your mouth
2. Press down the power button to switch on the device. Slightly close your mouth around the attachment to avoid the water splashing.
3. Position the water jet stream flow at a 90° angle to the teeth and gums (fig. 6 below).

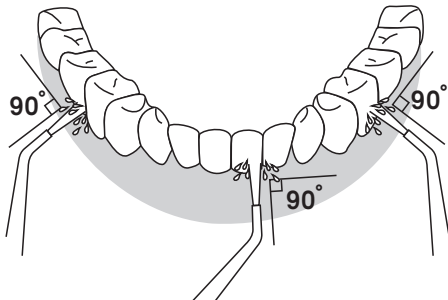4. By a short, repeated pressing of the button, the device switches between the Active, Sensitive and Impulse modes.
5. The unit remembers the last chosen mode and selects it the next time it is used.
6. The unit can be switched off by pressing the button at any time after 4 seconds of use.

7. Spot cleaning is used by pressing and holding the button. Point the attachment at the area to be cleaned. Press and hold the button when the unit is switched off. The device operates in sensitive mode. When the button is released, the unit switches off.

SPECIFIC FUNCTIONS ARE AS FOLLOWS:

Active mode:

High pressure water jet stream. The LED light glows once and stays illuminated during use.

Sensitive mode:

Gentle water jet stream. The LED light glows white twice and stays illuminated during use.

Pulse mode:

3 pulses per second. The LED light glows thrice and stays illuminated during use.

Spot Clean mode:

Continuous water jet stream while button is engaged. As soon as button is released, water stops.

HOW TO END USE

Auto shut-off turns the unit off after a 2-minute cycle is complete.

Turn it off by pressing the power button until the device is switched off.


RECOMMENDED TECHNIQUE

Cycle through the different modes by repeatedly pushing the power button.

The Flosserpik® Power Spa remembers the last chosen mode so that it does not have to be selected the next time you activate the device (**memory function**).

WATER JET


Slightly close lips to avoid splashing. Point the water jet stream at a 90° angle to your gum line.

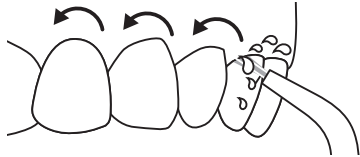Slowly move the attachment along the gum line. When finished, turn off the unit by pressing the power button on the base and empty the water tank.


ORTHODONTIC CLEANING

Using the ultrasoft brush tip, follow the recommended technique described below to clean around orthodontic brackets one tooth at a time.

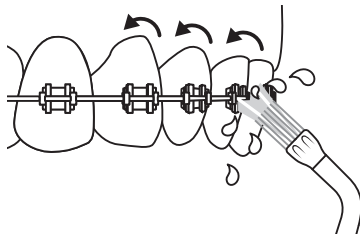BRUSHING WITH ULTRASOFT KONEX HD® TOOTHBRUSH


Insert brush head attachment into handle. Brush teeth as you normally would. After brushing, rinse mouth using the water jet stream.

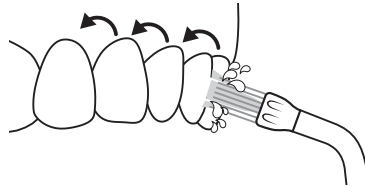SUBGINGIVAL ATTACHMENT

Designed for close comfort, the rubber tip delivers a concentrated jet stream of water or antibacterial solution that reaches deep into periodontal pockets.

Note: this tip is NOT intended for general use, and **MUST** be used with the Sensitive mode only. Place the soft tip against a tooth and gently slide the tip slightly under the gum line and into the pocket. Clean by gently following the gum line, inserting the tip between teeth. Repeat in all areas of the mouth.

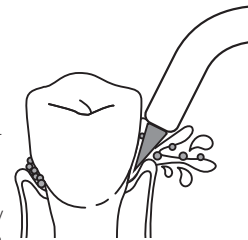RECOMMENDED REPLACEMENT INTERVALS FOR TIPS

Dentists recommend to replace the Brush tip every 3 months and the other tips every 6 months

USING MOUTHWASH AND OTHER SOLUTIONS

Your Flosserpiik® Power Spa can be used with non-alcoholic mouthwashes without essential oils. For general purpose, we recommend edel+white® Fresh+Protect mouthwash for natural freshness and antibacterial protection.

Important note: use of some solutions can shorten the life of your unit if it is not cleaned properly. You must thoroughly rinse out the unit after using any solution. Fill the tank with warm water to the lowest fill line, then aim the tip into the sink and run until the tank is empty.

CLEANING AND CARING FOR YOUR FLOSSERPIIK® POWER SPA ORAL WELLNESS SYSTEM

Cleaning the body:

Clean with a dry or damp towel.


Cleaning the water tank:

Wipe with a wet or dry towel and rinse clean with water if not using the device for a long period of time. Let it dry completely before storing.

Cleaning the attachments:

Use a wet or dry cloth to wipe or remove attachment, rinse with water and dry with a cloth before replacing.

Do not bend, pull or twist the suction pipe.

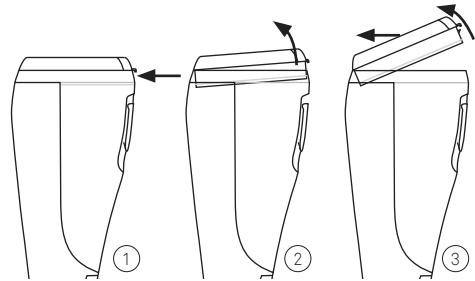 Open the hatch cover via the release button for cleaning (fig. 1-3). Never open when filling.

 Ensure hatch cover is properly sealed when reattaching (fig. 3-1).

 Gently wash the product with water before using. Clean the product every 6 months if it is not being used for a prolonged period. Clean the device with clean water or a little soap. If your water is hard, clean the appliance regularly with a little vinegar. Download instructions here (PDF): www.flosserpiik.com/care

 Do not use corrosive or abrasive detergents. Do not use hot water over 40°C. Do not put the water tank into dishwasher.

TROUBLESHOOTING GUIDE

The device doesn't start

If the device is newly purchased or has not been used for more than 3 months, charge for at least 4 hours. If the device still doesn't start, contact customer service.

The water pressure is low

The attachment is likely blocked. Replace attachment to increase water pressure. If this does not resolve the issue, contact customer service.

The attachment is deformed

If the attachment is damaged or defective in any way, replace it with a new one.

The battery can only be used for a short time after full charge

- The battery is at the end of its life. Contact customer service.
- The battery has not been sufficiently charged. Charge for at least 4 hours before use.

The charging indicator does not work

- The power is not properly connected. Ensure all cables and plugs are secure and switched on.
- There is a device circuit fault. Contact customer service.

There is no water coming out of the device

The device is newly purchased, unused for some time or there is a blockage in the attachment. Remove the attachment, flush the suction pipe and replace attachment.

TECHNICAL SPECIFICATIONS

This device complies with the guidelines of the EU (EMC/RoHS).

Rechargeable battery: Lithium-ion, ~30 uses in one full charge based on the 2-minute Active mode.

Charging time: ~6 hours

Charge mode: DC

Environment: 0-50°C (under standard atmospheric pressure), Applicable altitude below 5000m

The capacity of water tank: max.180ml

Water protection class: Unit: IPX7.

Material: ABS(bod), PC (water tank), PMMA(nozzle)

4x attachments

4 Modes: (Active/Sensitive/Pulse/Spot Clean)

WARRANTY

edel+white® guarantees its products for two years after the date of purchase. This warranty is the same in every country in which this device is sold by edel+white® or another authorised dealer. This warranty does not cover defects due to inappropriate use, regular wear or handling (especially of the interchangeable heads) and flaws that insignificantly affect value or functionality of the device. The warranty is invalidated if the device is repaired by unauthorised persons or if the body housing is opened.

In order to claim warranty, please send the entire device with proof of purchase to the official edel+white® partner in your country.

If you do not know your official edel+white® partner, please visit our homepage:

www.helpline.edelwhite.swiss

or contact us directly in Switzerland:

+41 (0)61 508 16 69

(The usual costs of a phone call into a foreign country will be charged).

DISPOSAL

The appliance must be disconnected from the power supply when removing the battery. Dispose of the discharged batteries correctly. Special containers are provided for disposing of rechargeable batteries at shops selling them or public collection points. When disposing the main unit, affix scotch tape or label to the electrodes of the battery. Do not throw the rechargeable battery in fire.

Correct disposal of this product:

This marking indicates that this product should not be disposed with other household wastes throughout the EU.

To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

Removal of battery and scrap processing - see end of manual

CUSTOMER SERVICE CENTRE

Unit: EW-WJ180 Flosserpik® Power Spa Art. 100-171

Replacement tips: EW-WJR2 Art. 100-172

For more information regarding this and other edel+white® products, visit

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



www.scanderra.com

SCRAP PROCESSING

Do not disassemble the product by yourself unless you want to dispose of it. Premature removal will damage the product. Operation is subject to basic safety precautions. Before removing the battery, you must unplug the power plug and charging plug.

1. Before removing the battery, make sure the battery is exhausted.
2. Remove the nozzle^A and the watertank^B in the direction of the arrow.
3. Remove the main body back cover in the direction of the arrow.
4. Use a screw driver to lift the 4 screws in the direction of the down arrow and remove the movement assembly.
5. Use a saw to cut off at the dotted line, and remove the battery. Property recycle or dispose of the battery. Be careful not to short the positive and negative terminals of the battery. Do not charge the product after removing the battery.

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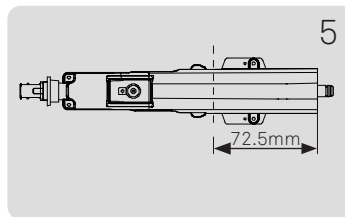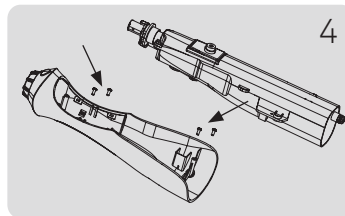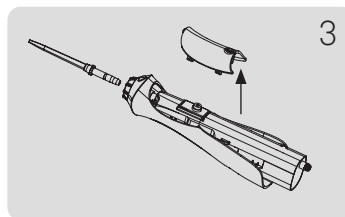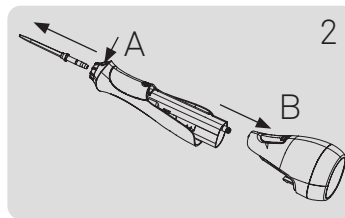EW-WJ180 Flosserpik® Power Spa | Art. 100-171

edel+white® and Flosserpik® are internationally registered trademarks of Scanderra GmbH, Basel, Switzerland. All rights and modifications reserved.